

Sensors and actuators used:

- Temperature: list in the tables from 5 to 7
- Heater: list in the tables from 5 to 7 and current leads heaters for magnet insert
- Level: LI670, LI680
- Valve: CV581

The user chooses:

- Temperature setpoint: list in the tables from 5 to 7
- Level: LI670mini, LI680mini
- Flow: FT581limit, FT583limit

Initial conditions:


- Sequences from 1 to 3 stopped
- Sequences 6 and 7 stopped
- Sequences from 10 to 20 stopped
- Sequence 8 in operation

This sequence drives the electrical heaters implicated in the warm-up of the cryostat. Each set of heaters has its own cycle. To limit the current draw when starting the heater, the sequences work this way: each set of heaters starts its cycle 10 seconds after the previous. The sequential starting of heaters may take several minutes. The cycles of heater control run in parallel until the user decides to stop the warm-up. The cycles then all stop at the same time.

The sequences used to warm-up the thermal shield and the helium circuits are very similar, but for the pressure vessel warm-up, the program checks the helium level. As long as the cryostat contains liquid helium, the heating is intermittent and allows evaporating the liquid stored in the cryostat. The heating operates in continuous mode when the cryostat is empty (Liquid helium level < 5%).

The sequence 22 can start only when the cryostat thermal shield cooling (sequence 6) and the cryostat helium cooling sequences are stopped.

The sequences 22.1 and 22.2 are started one after the other as described below:

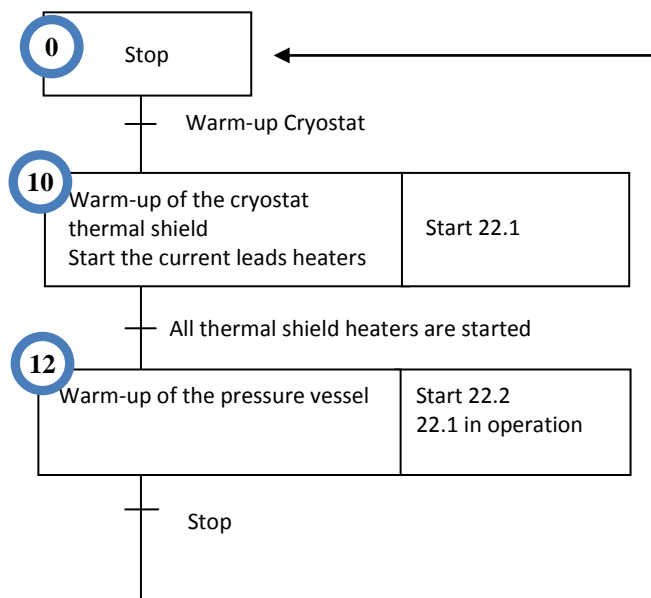

22-1: Warm up of the cryostat thermal shield

The sequence used for the warm-up of the cryostat thermal shield is similar to the valve box thermal shield sequence 21-1.

Vacuum or Liquid or Magnet mode

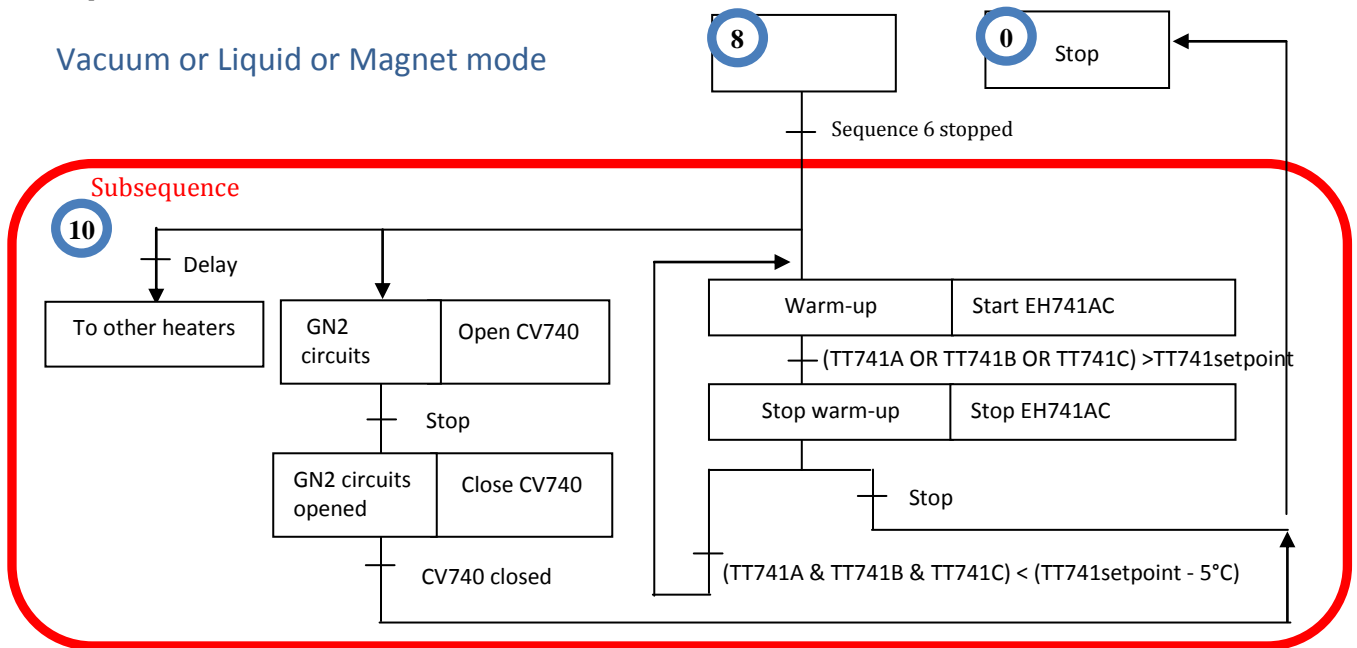


Table 5: Cryostat warm-up - Thermal shields

| Heater | Thermometer | Temperature setpoint |
|---------|-------------|----------------------|
| EH741AC | TT741AC | TT741setpoint |
| EH742AC | TT742AC | TT742setpoint |
| EH743AF | TT743AF | TT743setpoint |

22-2: Warm-up of the pressure vessel

This sequence 22-2 is similar to the sequence Warm up Valve Box 21-2.

Liquid mode

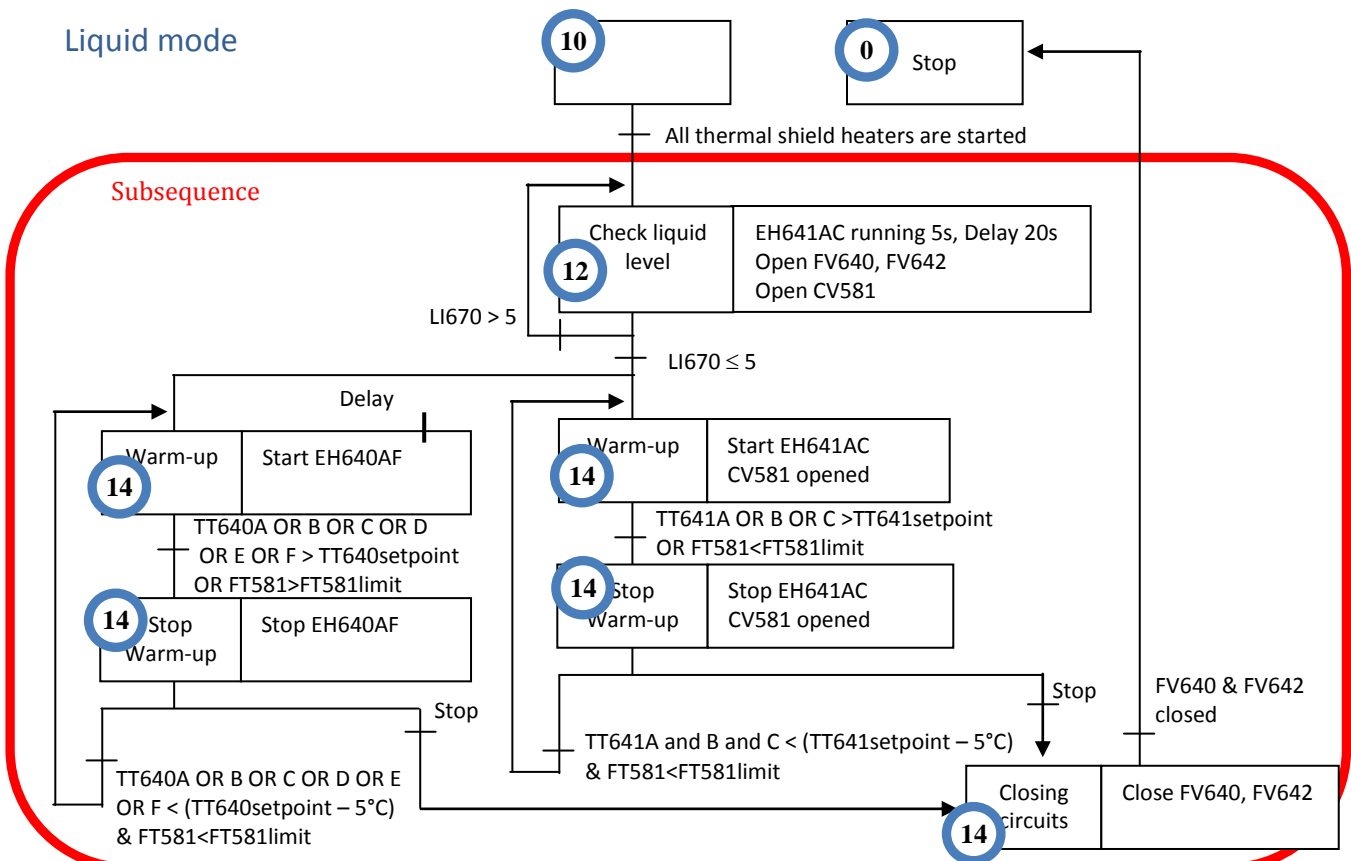


Table 6: Cryostat warm-up - Helium tank – Liquid mode

| Heater | Pt thermometer | Level | Temperature setpoint | CX thermometer | Temperature threshold |
|---------|----------------|-------|----------------------|----------------|-----------------------|
| EH640AF | TT640AF | LI670 | TT640setpoint | TT644 | 120K |
| EH641AC | TT641AC | LI670 | TT641setpoint | TT644 | 120K |

Magnet mode

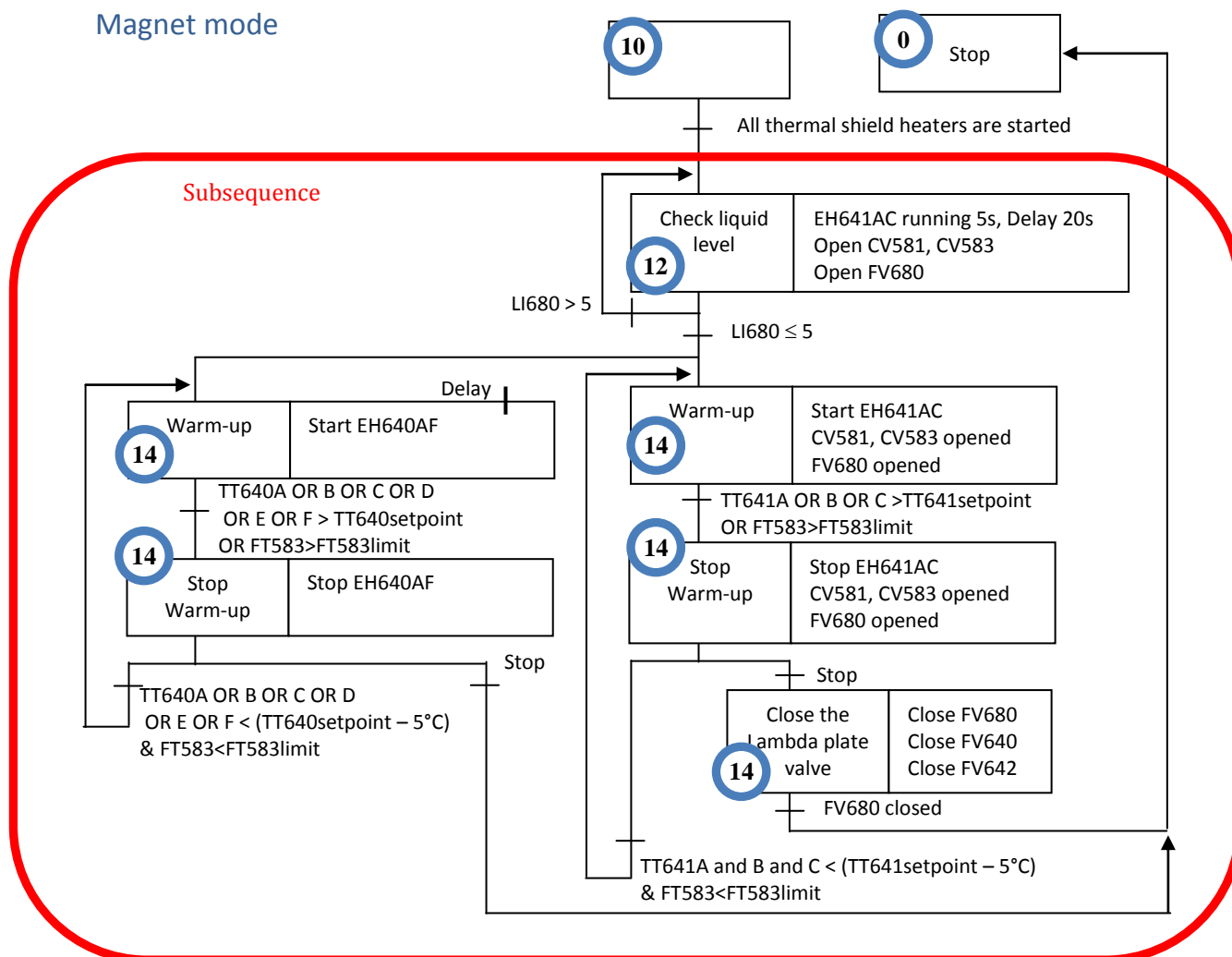


Table 7: Cryostat warm-up - Helium tank – Magnet mode

| Heater | Thermometer | Level | Temperature setpoint |
|---------|-------------|-------|----------------------|
| EH640AF | TT640AF | LI680 | TT640setpoint |
| EH641AC | TT641AC | LI680 | TT641setpoint |